

NHS Grampian Joint Wound Care Formulary

Product:

Flaminal® Forte

Product description:

Flaminal® is an enzyme alginogel®, used to debride, balance moisture levels and protect against bacteria as a biodegradable, soft and soothing wound interface reducing pain and trauma during wound care.

Indications for use:

Continuous debridement.
Kills bacteria without damaging the key healing cells.
Alginate balances moisture levels.
Protects wound edges.
Safe on all wound types at any healing stage.
Can be re-capped and re-used on the same patient.

Contraindications:

Should not be used where patients have had previous sensitivity to alginate dressings or Polyethylene glycol.
Should not be used in full thickness (third degree) burns.

Precautions:

None Known.
Flaminal® Forte cannot be used on the eyelids or in the eye. Should it come into contact with an eye, rinse the eye thoroughly with running water and consult a physician.

Product type:


Alginate gel containing antimicrobial enzymes, glucose oxidase and lactoperoxidase which exert their effect without damaging healthy skin cells. Debrides wound and manages moisture balance.

Product function:

Primary dressing
Moisture donating
Debriding

Wound type:

Can be used on a wide range of wounds at every stage of healing.

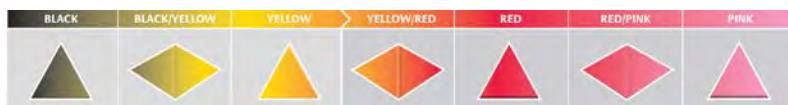

Product	Manufacturer	Available sizes	Pack size
Flaminal® Forte	Flen Pharma	15g	5
		50g	1